



CD8a Monoclonal Antibody (RPA-T8), eFluor™ 506, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Rat, Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), eFluor™ 506, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	RPA-T8
Conjugate	eFluor™ 506
Excitation/Emission Max	419/508 nm
Form	Liquid
Concentration	5 μL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2637468

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Flow Cytometry (Flow)	5 μL (0.125 μg)/test	24 Publications

Product Specific Information

Description: The RPA-T8 monoclonal antibody reacts with the human CD8a molecule, an approximately 32-34 kDa cell surface receptor expressed either as a heterodimer with the CD8 beta chain (CD8 alpha/beta) or as a homodimer (CD8 alpha/alpha). A majority of thymocytes and a subpopulation of mature T cells and NK cells express CD8a. CD8 binds to MHC class I and through its association with protein tyrosine kinase p56lck plays a role in T-cell development and activation of mature T cells.

Applications Reported: This RPA-T8 antibody has been reported for use in flow cytometric analysis.


Applications Tested: This RPA-T8 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 μ L (0.125 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


eFluor® 506 can be excited with the violet laser line (405 nm) and emits at 506 nm. We recommend using a 510/20 band pass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

Excitation: 405 nm; Emission: 506 nm; Laser: Violet Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD8a Monoclonal Antibody (RPA-T8), eFluor™ 506, eBioscience™

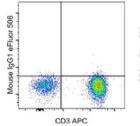

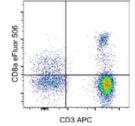

CD8a Antibody (69-0088-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD3 APC (Product # 17-0037-42) and 0.125 μg of Mouse IgG1 K Isotype Control eFluor® 506 (Product # 69-4714-82) (left) or 0.125 μg of Anti-Human CD8a eFluor® 506 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

□ 25 References

Immunohistochemistry (1)

Open life sciences

Changes of immune cells in patients with hepatocellular carcinoma treated by radiofrequency ablation and hepatectomy, a pilot study.

"Published figure using CD8a monoclonal antibody (Product # 69-0088-42) in Immunohistochemistry" Authors: Yu Z,Li G,Yu H,Asakawa T

Year 2022

Flow Cytometry (24)

Frontiers in immunology

Bispecific T-Cell Engaging Antibodies Against MUC16 Demonstrate Efficacy Against Ovarian Cancer in Monotherapy and in Combination With PD-1 and VEGF Inhibition.

"Published figure using CD8a monoclonal antibody (Product # 69-0088-42) in Flow Cytometry"

Authors: Yeku OO,Rao TD,Laster I,Kononenko A,Purdon TJ,Wang P,Cui Z,Liu H,Brentjens RJ,Spriggs D

Year 2021

Oncotarget

Multipeptide stimulated PBMCs generate T_{EM}/T_{CM} for adoptive cell therapy in multiple myeloma.

"Published figure using CD8a monoclonal antibody (Product # 69-0088-42) in Flow Cytometry" Authors: Vardam-Kaur T,Pathangey LB,McCormick DJ,Bergsagel PL,Cohen PA,Gendler SJ

Year 2021

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of substantiality for use in applications required by EPA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is immitted to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMINO PRODUCTS DURING HE WARRANTY PERFIOI DIS LIMITED. A REPLIAN FOR THE NON-CONFORMINO PRODUCTS, AT SELLER'S SOLE OPTION. THERE IS NO OBLICATION TO REPLIA, REPLACE OR REFUND FOR THE RON-CONFORMINO PRODUCTS, AT SELLER'S SOLE OPTION. THERE IS NO OBLICATION TO REPLIA, REPLACE OR REFUND FOR THE RON-CONFORMINO PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivo or in vivo therapeutic uses, or any type of