

CD134 (OX40) Monoclonal Antibody (OX-86), Super Bright 780, **eBioscience**

Product Details	
Size	100 μg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG1, kappa
Recommended Isotype Control	Rat IgG1 kappa Isotype Control (eBRG1), Super Bright 780, eBioscience
Class	Monoclonal
Туре	Antibody
Clone	OX-86
Conjugate	Super Bright 780
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage Conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2734876

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	1.0 µg/test	-

Product Specific Information

Description: The OX-86 monoclonal antibody reacts with mouse CD134, also known as OX40. A member of the TNF receptor superfamily, CD134 is a 50 kDa type I membrane glycoprotein expressed by activated mouse T lymphocytes. Rat CD134 was initially identified as an activation marker only on activated CD4+ T cells. In contrast, mouse CD134 is expressed by both activated CD4+ and CD8+ T cells. The interaction of CD134 with CD252 (OX40 ligand) has been implicated in T cell-dependent humoral responses, regulation of primary T cell expansion, survival of T cells, size of the memory T cell pool, and regulation of tolerance in the CD4+ T cell compartment.

Applications Reported: The OX-86 antibody has been reported for use in flow cytometric analysis.

Applications Tested: The OX-86 antibody has been tested by flow cytometric analysis of activated mouse splenocytes. This can be used at less than or equal to 1.0 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 μL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Super Bright 780 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 780 nm. We recommend using a 780/60 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

In some experiments, we have observed that compensation values for Super Bright 780-conjugated antibodies are higher in the violet 450/50 channel when using UltraComp eBeads microspheres (Product # 01-2222-42) as compared to single-color stained cells. In such circumstances, we would recommend setting compensation with cells. We have also observed this in some experiments using AbC Total Antibody Compensation beads (Product # A10497).


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 405 nm; Emission: 780 nm; Laser: Violet Laser.

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD134 (OX40) Monoclonal Antibody (OX-86), Super Bright 780, eBioscience

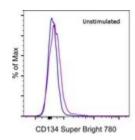

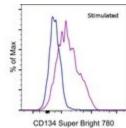

CD134 (OX40) Antibody (78-1341-82) in Flow

Staining of unstimulated (left) or 3-day Con A (00-4978-93)-stimulated C57Bl/6 mouse splenocytes (right) with 1.0 μ g of Rat IgG1 kappa Isotype Control, Super Bright 780 (Product # 78-4301-82) (blue histogram) or 1.0 μ g of CD134 (OX40) Monoclonal Antibody, Super Bright 780 (purple histogram). Total viable cells, as determined by Fixable Viability Dye eFluor® 450 (Product # 65-0863-14), were used for analysis.

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYERS EXCLUSIVE REMEDY FOR NON-CORNING PRODUCTS ON THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OR REPLIAND FOR THE NON-CONFORMING PRODUCTS, AT SELECES'S SOLE OPTION. THERE IS NO. OD OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption to human or animals.