



IL-8 (CXCL8) Polyclonal Antibody


Product Details	
Size	100 μg
Species Reactivity	Human
Published Species	Human
Host/Isotype	Rabbit / IgG
Class	Polyclonal
Туре	Antibody
Conjugate	Unconjugated
Immunogen	Highly purified (>98% pure) recombinant human endothelial IL-8 (77 amino acid residues).
Form	Lyophilized
Purification	purified
Storage buffer	PBS, pH 7.2
Contains	no preservative
Storage conditions	-20°C
RRID	AB_2536263

Applications	Tested Dilution	Publications
Western Blot (WB)	1:1,000	-
Immunohistochemistry (IHC)	Assay-dependent	8 Publications
Immunohistochemistry (Paraffin) (IHC (P))	-	5 Publications
Immunohistochemistry (Frozen) (IHC (F))	-	2 Publications
Immunocytochemistry (ICC/IF)	-	1 Publication
ELISA (ELISA)	Assay-dependent	1 Publication
Neutralization (Neu)	Assay-dependent	1 Publication
In vitro Assay (IV)	-	1 Publication

Product Specific Information

Centrifuge vial prior to opening. Reconstitute in sterile water to 0.1-1.0 mg/mL.

Product Images For IL-8 (CXCL8) Polyclonal Antibody

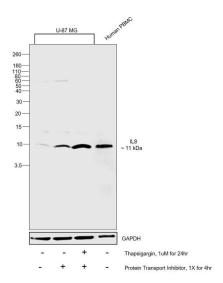

IL-8 (CXCL8) Antibody (AHC0881) in WB

Western blot was performed using Anti-IL-8 (CXCL8) Polyclonal Antibody (Product # AHC0881) and a 11kDa band corresponding to IL-8 (CXCL8) was observed in Human PBMC and U-87 MG with Thapsigargin and PTI treatment. Whole cell extracts (30 µg lysate) of U-87 MG (Lane 1), U-87 MG treated with PTI (1X for 4hr) (Lane 2), U-87 MG treated with Thapsigargin (1uM for 24hr) followed by PTI (1X for 4hr) (Lane 3) and Human PBMC (Lane 4) were electrophoresed using NovexTM 16% Tricine Protein Gel (Product # EC6695BOX). Resolved proteins were then transferred onto a Nitrocellulose membrane (Product # IB23001) by iBlot® 2 Dry Blotting System (Product # IB21001). The blot was probed with the primary antibody (1:1000 dilution) and detected by chemiluminescence with Goat anti-Rabbit IgG (Heavy Chain) SuperclonalTM Recombinant Secondary Antibody, HRP (Product # A27036, 1: 4000 dilution) using the iBright FL 1000 (Product # A32752). Chemiluminescent detection was performed using Novex® ECL Chemiluminescent Substrate Reagent Kit (Product # WP20005).

View more figures on thermofisher.com

□ 19 References

Immunohistochemistry (8)

Clinical and experimental gastroenterology

Toll-Like Receptor 3 as a Recurrence Risk Factor and a Potential Molecular Therapeutic Target in Colorectal Cancer.

"AHC0881 was used in Immunohistochemistry to determine that TLR3 expression can be used as a risk factor for the recurrence of colorectal cancer, as well as a potential target for new therapeutic strategies."

Authors: Yoshida T,Miura T,Matsumiya T,Yoshida H,Morohashi H,Sakamoto Y,Kurose A,Imaizumi T,Hakamada K

Year 2022

Species Human

Dilution 1:100

Diagnostics (Basel, Switzerland)

Prognostic Significance of Lymphocyte Infiltration and a Stromal Immunostaining of a Bladder Cancer Associated Diagnostic Panel in Urothelial Carcinoma.

"AHC0881 was used in Immunohistochemistry to report on the presence of 10 biomarkers in the stroma of human bladder cancer and to determine the clinical significance of the signature."

Authors: Furuya H,Chan OTM,Hokutan K,Tsukikawa Y,Chee K,Kozai L,Chan KS,Dai Y,Wong RS,Rosser CJ

Year 2019

Species Human

Dilution 1:200

View more IHC references on thermofisher.com

Immunohistochemistry (Paraffin) (5)

BioMed research international

Immunohistochemical analysis of IL-6, IL-8/CXCR2 axis, Tyr p-STAT-3, and SOCS-3 in lymph nodes from patients with chronic lymphocytic leukemia: correlation between microvascular characteristics and prognostic significance.

"AHC0881 was used in immunohistochemistry - paraffin section to study the correlation of microvascular characteristics and prognostic differences with immunohistochemical analysis of IL-6, IL-8/CXCR2 axis, Tyr p-STAT-3, and SOCS-3 in chronic lymphocytic le"

Authors: Levidou G,Sachanas S,Pangalis GA,Kalpadakis C,Yiakoumis X,Moschogiannis M,Sepsa A,Lakiotaki E,Milionis V,Kyrtsonis MC,Vassilakopoulos TP,Tsirkinidis P,Kontopidou F,Kokoris S,Siakantaris M,Angelopoulou M,Papadaki H, Kavantzas N,Panayiotidis P,Patsouris E,Korkolopoulou P

Year 2015

Diagnostic pathology

Validation and clinicopathologic associations of a urine-based bladder cancer biomarker signature.

"AHC0881 was used in immunohistochemistry - paraffin section to validate the expression of a urine-based bladder cancer associated diagnostic signature"

Authors: Zhang G,Gomes-Giacoia E,Dai Y,Lawton A,Miyake M,Furuya H,Goodison S,Rosser CJ

Year 2014

Species Human

Dilution 1:200

View more IHC (P) references on thermofisher.com

More applications with references on thermofisher.com

IHC (F) (2) ICC/IF (1) ELISA (1) Neu (1) IV (1)

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and carcompanying package insenser (Documentation). No claim of suitability for use in applications required by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is immitted to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE EMBLED'FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERFIOID IS LIMITED. A REFUND FOR THE NON-CONFORMING PRODUCTS (S) AT SELLER'S SOLE OPTION. THERE IS NO AGAIN TO TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, SOLE OPTION. THERE IS NO AGAIN TO THE REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS (S) AT SELLER'S SOLE OPTION. THERE IS NO AGAIN TO THE REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLICENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized c