CD43 Monoclonal Antibody (eBio84-3C1 (84-3C1)), NovaFluor™ Yellow 570, eBioscience™

Product Details

Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Class	Monoclonal
Туре	Antibody
Clone	eBio84-3C1 (84-3C1)
Conjugate	NovaFluor™ Yellow 570
Excitation/Emission Max	551/566 nm
Form	Liquid
Concentration	0.2 µg/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.2 μg)/test	-

Product Specific Information

Description: The monoclonal antibody ebio84-3C1 reacts with CD43, which is also known as leukosialin, galactoglycoprotein and sialophorin. CD43 is a sialomucin which like many mucins can have both adhesive and anti-adhesive functions. Expression of CD43 is found on most leukocytes except resting B lymphocytes. Proteolytic processing upon activation decreases surface expression. CD43 is involved in activation of T cells, B cells, NK cells, and monocytes. The counter-receptor for CD43 is CD169/SIGLEC-1, which is expressed on macrophages.

Each product contains 1 vial of NovaFluor conjugate and 1 vial of CellBlox Plus Blocking Buffer .

Applications Reported: This ebio84-3C1 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This ebio84-3C1 antibody has been pre-diluted and tested by flow cytometric analysis of normal human peripheral blood cells. This may be used at 5 μ L (0.2 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test.

NovaFluor dyes are not compatible with DNA intercalating viability dyes. Do not use viability dyes such as propidium iodide, 7actinomycin D (7-AAD) and DAPI. Invitrogen LIVE/DEAD Fixable Dead Cell stains are recommended for use with NovaFluor dyes.


This NovaFluor conjugate has been updated to ship with CellBlox Plus Blocking Buffer (Cat. No. (C001T06F01)). This buffer contains formulation improvements over CellBlox. CellBlox Plus Blocking Buffer is required for optimal staining with NovaFluor conjugates and should be used in all experiments where NovaFluor conjugates are used. Whenever possible, we recommend adding CellBlox Plus Blocking Buffer to antibody cocktails/master mixes prior to combining with cells. Add 5 μ L per sample (regardless of the number of NovaFluors in your panel) to use the antibody cocktail as intended. For single-color controls, use 5 μ L of CellBlox Blocking Buffer per 100 μ L of cell sample containing 10^3 to 10^8 cells.

1

NovaFluor conjugates are based on Phiton[™] technology utilizing novel nucleic acid dye structures that allow for engineered fluorescent signatures with consideration for spillover and spread impacts. Learn more


Excitation: 568 nm; Emission: 568 nm; Laser: 561 nm (Yellow) Laser

Product Images For CD43 Monoclonal Antibody (eBio84-3C1 (84-3C1)), NovaFluor™ Yellow 570, eBioscience™

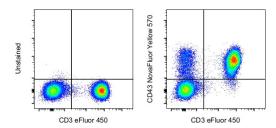


CD43 Antibody (H069T03Y01-A) in Flow

Normal human peripheral blood cells were unstained (left) or stained with CD43 Monoclonal Antibody, NovaFluor Yellow 570 (right). All cells were co-stained with CD3 Monoclonal Antibody, eFluor 450 (Product # 48-0038-82). Total viable cells in the lymphocyte gate were used for analysis, as determined by LIVE/DEAD Blue (Product # L34962). Data was acquired on a 5-laser Cytek Aurora and unmixed with autofluorescence extraction.

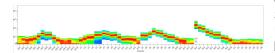

CD43 Antibody (H069T03Y01-A) in Flow

Spectral signature for NovaFluor Yellow 570 collected on a 5-laser Cytek Aurora Full Spectrum flow cytometer using Cytek assay settings. Human peripheral blood mononuclear cells were stained with anti-human CD4 (SK3) and signatures displayed following gating on the lymphocyte population.

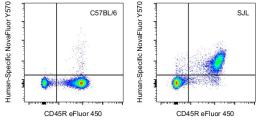

CD43 Antibody (H069T03Y01-A) in Flow

NovaFluor Yellow 570 non-specific staining of B cells in the SJL strain of mice. Splenocytes from C57BL/6 (left) and SJL (right) strains of mice were stained with Anti-Mouse CD45R (B220) Monoclonal Antibody conjugated to eFluor 450 and a non-cross-reactive, human-specific monoclonal antibody conjugated to NovaFluor Yellow 570. These data show that NovaFluor Yellow 570-conjugated antibodies non-specifically stain B cells in SJL mice (right). Non-specific staining has not been observed in C57BL/6 (left) mice and BALB/c mice (data not shown).

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and guarding package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is imited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES SOF ANY PARTICLAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE EMEMORY FOR NON-CONFERNING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO OBLICATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO OBLICATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO OBLICATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION TO THER WARRANT OF FOR CHAILED THE WARRANT OF FOR CHAILED THE IN OUR CHAILED THE WARRANT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS. UNess otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including w