Performance guarenteed'

# Clathrin Heavy Chain Monoclonal Antibody (X22), Alexa Fluor™ 488

#### **Product Details**

Size	50 µL
Species Reactivity	Bovine, Hamster, Human, Mouse, Non-human primate, Rat
Host/Isotype	Mouse / IgG1
Class	Monoclonal
Туре	Antibody
Clone	X22
Conjugate	Alexa Fluor™ 488
Excitation/Emission Max	499/520 nm
Immunogen	Purified human brain clathrin heavy chain
Form	Liquid
Concentration	1 mg/mL
Purification	Protein A
Storage buffer	PBS, pH 7.4, with 1mg/mL BSA
Contains	0.05% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2662916

Applications	Tested Dilution	Publications
Western Blot (WB)	-	4 Publications
Immunocytochemistry (ICC/IF)	1-10 μg/mL	3 Publications
Immunoprecipitation (IP)	-	2 Publications


#### **Product Specific Information**

MA1065A488 detects clathrin heavy chain in non-human primate, bovine, human, rat and mouse tissues as well as hamster (CHO) cells.

The MA1065A488 immunogen is purified human brain clathrin heavy chain. Electron microscopy and proteolysis mapping demonstrate that MA1065A488 binding occurs towards the central hub of the triskelion, N-terminal to the light chain binding regions.

1

### Product Images For Clathrin Heavy Chain Monoclonal Antibody (X22), Alexa Fluor™ 488

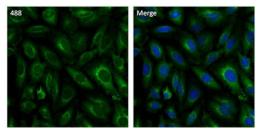

#### Clathrin Heavy Chain Antibody (MA1-065-A488) in ICC/IF

Immunofluorescent analysis of Clathrin Heavy Chain (green) in HeLa cells. The cells were fixed with 4% Paraformaldehyde in PBS for 15 minutes at room temperature, and blocked with 3% BSA in PBS (Product # 37525) for 30 minutes at room temperature. Cells were stained with a Clathrin Heavy Chain Monoclonal Antibody, AlexaFluor 488 conjugate (Product # MA1-065-A488) at a dilution of 5 µg/mL in blocking buffer for 1 hour at room temperature protected from light. Nuclei (blue) were stained with Hoechst Dye (Product # 62249) at a dilution of 1: 10,000 in blocking buffer. Images were taken on a Thermo Scientific ToxInsight Instrument at 20X magnification.

View more figures on thermofisher.com

#### 9 References

#### Western Blot (4)

cells induces inflammatory cyte	ression in metastatic intestinal epithelial kine production and promotes
invasiveness.	
"Published figure using Clathrin Heavy Chain mon	clonal antibody (Product # MA1-065-A488) in Western Blot"
Authors: Bugge M,Bergstrom B,Eide OK,Solli H,Kj	anstad IF,Stenvik J,Espevik T,Nilsen NJ
Scientific reports	
POCK1 is a noval Pact offector	to regulate tubular endocytic membrane
formation during clathrin-indep	to regulate tubular endocytic membrane

"Published figure using Clathrin Heavy Chain monoclonal antibody (Product # MA1-065-A488) in Western Blot" Authors: Soriano-Castell D,Chavero A,Rentero C,Bosch M,Vidal-Quadras M,Pol A,Enrich C,Tebar F

#### Immunocytochemistry (3)

View more WB references on thermofisher.com

**Year** 2017

Year 2017

Frontiers in molecular neuroscience	<b>Year</b> 2023		
An Aberrant Phosphorylation of Amyloid Precursor Protein Tyrosine			
Regulates Its Trafficking and the Binding to the Clathrin Endocytic			
Complex in Neural Stem Cells of Alzheimer's Disease Patients.			
Published figure using Clathrin Heavy Chain monoclonal antibody (Product # MA1-065-A488) in Immunoprecipitation"			
Authors: Poulsen ET,Iannuzzi F,Rasmussen HF,Maier TJ,Enghild JJ,Jørgensen AL,Matrone C			
Nature neuroscience	Year		
Control role for DICALM in emulaid, blood brein berrier treneouteric and			

# Central role for PICALM in amyloid- blood-brain barrier transcytosis and clearance.

"Published figure using Clathrin Heavy Chain monoclonal antibody (Product # MA1-065-A488) in Immunofluorescence" Authors: Zhao Z,Sagare AP,Ma Q,Halliday MR,Kong P,Kisler K,Winkler EA,Ramanathan A,Kanekiyo T,Bu G,Owens NC,Rege SV,Si G,Ahuja A,Zhu D,Miller CA,Schneider JA,Maeda M,Maeda T,Sugawara T,Ichida JK,Zlokovic BV

View more ICC/IF references on thermofisher.com

## More applications with references on thermofisher.com

# IP (2)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set for hin the Production documentation, specifications and/or cascompanying package inserts ("Documentation in contrast, provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is enclined to one year from date of shipment when the Product is subjected to normal, proper and individuals. Unless otherwise stated in the Documentation, this warranty is one system that any Product will conform to such model or sample (unnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INPELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTBALLTY, FITNESS FOR ANY PARTICUAR PURPOSE, OR NON INFRINCEMENT. BUYER'S SCILUSVIE REMEUPS FOR NON-CONFERNING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BULGATION TO REPARIN, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product in the documentation and rainals.