



## **FOXP3 Polyclonal Antibody**

<b>Product Details</b>	
Size	100 μL
Species Reactivity	Bovine, Human, Mouse, Rat
Host/Isotype	Rabbit / IgG
Class	Polyclonal
Туре	Antibody
Conjugate	Unconjugated
Immunogen	A peptide containing amino acids 43-100 of mouse FOXP3
Form	Liquid
Concentration	1 mg/mL
Purification	Antigen affinity chromatography
Storage buffer	PBS
Contains	0.05% sodium azide
Storage conditions	-20°C or -80°C if preferred
RRID	AB_2278500

Applications	Tested Dilution	Publications
Western Blot (WB)	1:2,000	-
Immunohistochemistry (IHC)	-	2 Publications
Immunohistochemistry (Paraffin) (IHC (P))	1:400	-
Immunocytochemistry (ICC/IF)	1:500	-
Flow Cytometry (Flow)	1:8,000	-
Immunoprecipitation (IP)	Assay-Dependent	-


## **Product Specific Information**

Prior to immunostaining paraffin tissues, antigen retrieval with sodium citrate buffer (pH 6.0) is recommended.

Suggested positive control: mouse CD4+CD25+ T cells.


Epitope affinity purified.

## **Product Images For FOXP3 Polyclonal Antibody**

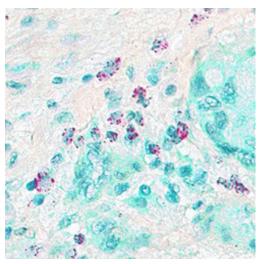


#### FOXP3 Antibody (PA1-46126) in IHC (P)

Immunohistochemical analysis of FOXP3 in formalin-fixed paraffin-embedded tissue sections of human breast cancer. Samples were incubated in FOXP3 polyclonal antibody (Product # PA1-46126) using a dilution of 1:400. ACDs Integrated Co-Detection Workflow was performed using ACD RNAScope Probe. Tissue was stained on Leica Bond RX using RNAscope (TM) 2.5 LS Reagent Kit-RED, BOND Polymer Refine Detection (DAB) and Hematoxylin, BOND Polymer Refine Red Detection and Hematoxylin and RNAscope (TM) 2.5 LS Green Accessory Pack. Tissue was counterstained with 50% hematoxylin (blue). CD8A mRNA (red) and FOXP3 protein (green).

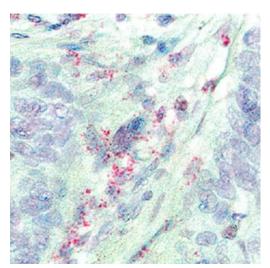

## FOXP3 Antibody (PA1-46126) in IHC (P)

Immunohistochemical analysis of FOXP3 in formalin-fixed paraffin-embedded tissue sections of human breast cancer. Samples were incubated in FOXP3 polyclonal antibody (Product # PA1-46126) using a dilution of 1:200. ACDs Integrated Co-Detection Workflow was performed using ACD RNAScope Probe. Tissue was stained on Leica Bond RX using RNAscope (TM) 2.5 LS Reagent Kit-RED, BOND Polymer Refine Detection (DAB) and Hematoxylin, BOND Polymer Refine Red Detection and Hematoxylin and RNAscope (TM) 2.5 LS Green Accessory Pack. Tissue was counterstained with 50% hematoxylin (blue). CD4 mRNA (red) and FOXP3 protein (green).

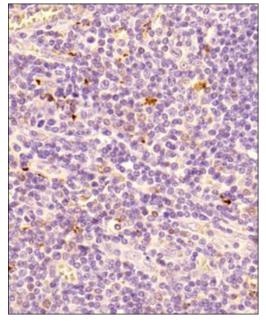

## FOXP3 Antibody (PA1-46126) in IHC (P)

Immunohistochemical analysis of FOXP3 in Tissue section of human tonsil. Samples were incubated in FOXP3 polyclonal antibody (Product # PA1-46126) using a dilution of 1:200 followed by HRP-labeled anti-rabbit IgG secondary antibody and DAB reagent. Nuclei of cells were counterstained with hematoxylin. This FoxP3 antibody generated an expected and a specific nuclear staining of FOXP3 in a subset of cells which are potentially the Treg cells.

View more figures on thermofisher.com

#### **□ 2 References**

### Immunohistochemistry (2)

#### Oncoimmunology

## Comedications influence immune infiltration and pathological response to neoadjuvant chemotherapy in breast cancer.

"Published figure using FOXP3 polyclonal antibody (Product # PA1-46126) in Immunofluorescence"

Authors: Hamy AS,Derosa L,Valdelièvre C,Yonekura S,Opolon P,Priour M,Guerin J,Pierga JY,Asselain B,De Croze D, Pinheiro A,Lae M,Talagrand LS,Laas E,Darrigues L,Grandal B,Marangoni E,Montaudon E,Kroemer G,Zitvogel L,Reyal F

**Year** 2021

#### Cancer medicine

# p65/miR-23a/CCL22 axis regulated regulatory T cells recruitment in hepatitis B virus positive hepatocellular carcinoma.

"Published figure using FOXP3 polyclonal antibody (Product # PA1-46126) in Immunofluorescence" Authors: Li ZQ,Wang HY,Zeng QL,Yan JY,Hu YS,Li H,Yu ZJ

**Year** 2020

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CORNORNING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED AND PREVENT OF OR PREVENT FOR THE NON-CONFORNING PRODUCTS; AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORNING PRODUCTS; AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORNING PRODUCTS; AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORNING PRODUCTS. AT SHAPP AND PRODUCTS OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivo or in vivo therapeutic use