

Calcium Sensing Receptor Polyclonal Antibody

Product Details	
Size	100 µL
Species Reactivity	Human, Rat
Published Species	Human, Mouse
Host/Isotype	Rabbit / IgG
Class	Polyclonal
Type	Antibody
Conjugate	Unconjugated
Immunogen	Synthetic peptide corresponding to residues A(12) L A W H S S A Y G P D Q R A Q(27) of rat CaSR.
Form	Liquid
Concentration	1 mg/mL
Purification	Antigen affinity chromatography
Storage buffer	PBS
Contains	0.05% sodium azide
Storage conditions	-20° C, Avoid Freeze/Thaw Cycles
RRID	AB_2071477


Applications	Tested Dilution	Publications
Western Blot (WB)	1:100	3 Publications
Immunohistochemistry (IHC)	-	6 Publications
Immunohistochemistry (Paraffin) (IHC (P))	-	1 Publication
Immunocytochemistry (ICC/IF)	1:10-1:100	2 Publications

Product Specific Information

PA1-934A detects calcium sensing receptor (CaSR) from human and rat tissues.

PA1-934A has been successfully used in Western blot and ICC/IF procedures. By Western blot, this antibody detects an ~120 kDa protein representing CaSR from HEK393 cells and rat brain homogenate.


The PA1-934A immunizing peptide corresponds to amino acid residues 12-27 from rat CaSR. This peptide (Cat.# PEP-009) is available for use in neutralization and control experiments.

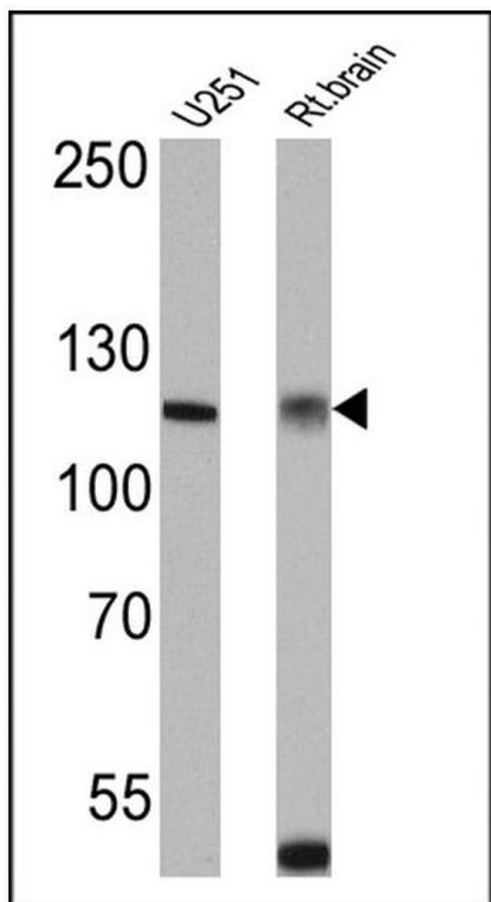


Calcium Sensing Receptor Antibody (PA1-934A) in WB

Western blot analysis of Calcium Sensing Receptor was performed by loading 25 µg of U251 (lane 1) and rat brain (lane 2) cell lysates onto an SDS polyacrylamide gel. Proteins were transferred to a PVDF membrane and blocked at 4°C overnight. The membrane was probed with a Calcium Sensing Receptor polyclonal antibody (Product # PA1-934A) at a dilution of 1:500 overnight at 4°C, washed in TBST, and probed with an HRP-conjugated secondary antibody for 1 hr at room temperature in the dark. Chemiluminescent detection was performed using Pierce ECL Plus Western Blotting Substrate (Product # 32132). Results show a band at ~120 kDa.

Calcium Sensing Receptor Antibody (PA1-934A) in ICC/IF

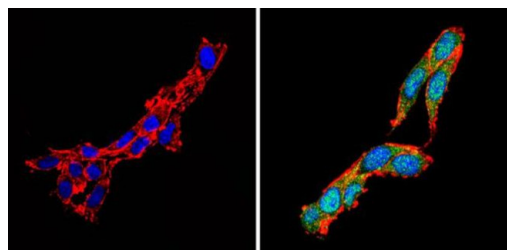

Immunofluorescent analysis of Calcium Sensing Receptor (green) showing staining in the cytoplasm and nucleus of PC12 cells (right) compared to a negative control without primary antibody (left). Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with a Calcium Sensing Receptor polyclonal antibody (Product # PA1-934A) in 3% BSA-PBS at a dilution of 1:20 and incubated overnight at 4 °C in a humidified chamber. Cells were washed with PBST and incubated with a DyLight-conjugated secondary antibody in PBS at room temperature in the dark. F-actin (red) was stained with a fluorescent red phalloidin and nuclei (blue) were stained with Hoechst or DAPI. Images were taken at a magnification of 60x.

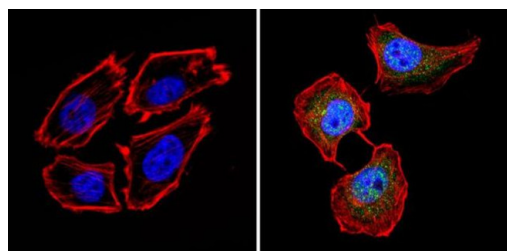

Calcium Sensing Receptor Antibody (PA1-934A) in ICC/IF

Immunofluorescent analysis of Calcium Sensing Receptor (green) showing staining in the cytoplasm and nucleus of U251 cells (right) compared to a negative control without primary antibody (left). Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with a Calcium Sensing Receptor polyclonal antibody (Product # PA1-934A) in 3% BSA-PBS at a dilution of 1:20 and incubated overnight at 4 °C in a humidified chamber. Cells were washed with PBST and incubated with a DyLight-conjugated secondary antibody in PBS at room temperature in the dark. F-actin (red) was stained with a fluorescent red phalloidin and nuclei (blue) were stained with Hoechst or DAPI. Images were taken at a magnification of 60x.

View more figures on thermofisher.com

Western Blot (3)

Journal of cell science	Year 2019
Evidence for a regulated Ca²⁺ entry in proximal tubular cells and its implication in calcium stone formation.	
"PA1-934A was used in Immunohistochemistry to show the existence of a regulated transcellular Ca2+ entry pathway in luminal membrane proximal tubule cells induced by Ca2+ sensing receptor-mediated activation of transient receptor potential canonical 3 channels."	
Authors: Ibeh CL,Yiu AJ,Kanaras YL,Paal E,Birnbaumer L,Jose PA,Bandyopadhyay BC	
European journal of endocrinology	Year 2012
A novel germline inactivating mutation in the CASR gene in an Italian kindred affected by familial hypocalciuric hypercalcemia.	Species Human
"PA1-934A was used in western blot to study a novel germline CASR mutation in a case of familial hypocaliuric hypercalcemia in an Italian family"	Dilution 1:500
Authors: Falchetti A,Gozzini A,Terranegra A,Soldati L,Vezzoli G,Leoncini G,Giusti F,Franceschelli F,Masi L,Tanini A,Cavalli L,Brandi ML	

View more WB references on thermofisher.com

Immunohistochemistry (6)

Journal of cell science	Year 2019
Evidence for a regulated Ca²⁺ entry in proximal tubular cells and its implication in calcium stone formation.	
"PA1-934A was used in Immunohistochemistry to show the existence of a regulated transcellular Ca2+ entry pathway in luminal membrane proximal tubule cells induced by Ca2+ sensing receptor-mediated activation of transient receptor potential canonical 3 channels."	
Authors: Ibeh CL,Yiu AJ,Kanaras YL,Paal E,Birnbaumer L,Jose PA,Bandyopadhyay BC	
Transplantation proceedings	Year 2013
Decreased parathyroid Klotho expression is associated with persistent hyperparathyroidism after kidney transplantation.	Species Human
"PA1-934A was used in immunohistochemistry to study persistent hypoparathyroidism in kidney transplant patients and the role of reduced parathyroid expression of Klotho"	
Authors: Hong YA,Choi DE,Lim SW,Yang CW,Chang YK	

View more IHC references on thermofisher.com

More applications with references on thermofisher.com

- IHC (P) (1)
- ICC/IF (2)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.