

NOV Polyclonal Antibody

Product Details

Size	100 μL	
Species Reactivity	Human	
Host/Isotype	Rabbit / IgG	
Class	Polyclonal	
Туре	Antibody	
Conjugate	Unconjugated	
Immunogen	Recombinant protein corresponding to Human NOV	
Form	Liquid	
Concentration	0.3 mg/mL	
Purification	Antigen affinity chromatography	
Storage buffer	PBS, pH 7.2, with 40% glycerol	
Contains	0.02% sodium azide	
Storage conditions	Store at 4°C short term. For long term storage, store at -20°C, avoiding freeze/thaw cycles.	
RRID	AB_2644755	

Applications	Tested Dilution	Publications
Western Blot (WB)	0.04-0.4 µg/mL	-
Immunohistochemistry (IHC)	1:200-1:500	-
Immunocytochemistry (ICC/IF)	0.25-2 µg/mL	-


Product Specific Information

Immunogen sequence: SCCLVCARQR GESCSDLEPC DESSGLYCDR SADPSNQTGI CTAVEGDNCV FDGVIYRSGE KFQPSCKFQC TCRDGQIGC

Highest antigen sequence identity to the following orthologs: Mouse - 81%, Rat - 80%.

1

Product Images For NOV Polyclonal Antibody

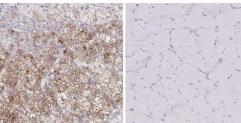

NOV in Adrenal gland


NOV in Skeletal muscle

RNA [TPM]

RNA [TPM] 200 400 600 800 1000


NOV Antibody (PA5-54055)

Relative expression in different tissues in IHC: Detection of differential expression levels of NOV demonstrates antibody specificity. Immunohistochemical analysis of NOV using anti-NOV Polyclonal Antibody (Product # PA5-54055), shows significant staining of NOV in human adrenal gland and shows minimal or weak staining in human skeletal muscle tissues. The relative expression levels of NOV within each tissue is shown using RNA-Seq. {RE}

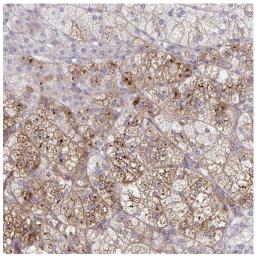

NOV Antibody (PA5-54055) in IHC

Immunohistochemical staining of NOV in human adrenal gland using NOV Polyclonal Antibody (Product # PA5-54055) shows moderate to strong cytoplasmic positivity in glandular cells.

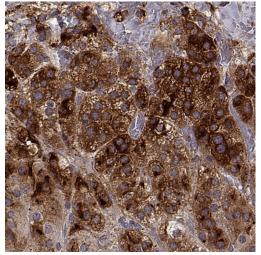

NOV Antibody (PA5-54055) in IHC

Immunohistochemical staining of NOV in human adrenal gland using a NOV Polyclonal Antibody (Product # PA5-54055) shows strong cytoplasmic positivity in zona glomerulosa.

View more figures on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty are provided herein is valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty are provided herein is valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty are provided herein is valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLED WARRANTES OF MERCHANTABLITY, FITHESS FOR ANY PARTICULA PURPOSE, OR NON INFRINCEMENT. BUYER'S SCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACE OR OR REFUND FOR THE NOD-CONFORMING PRODUCTS ON DANADER FOR WITHOUT TO FOR CE MALEVINE, (I) MISUSE, FAULT OR NEGLEGANCE OR BY BUYER, (II) USE OF THE PRODUCTS IN A MANNER FOR WITHOUT THEY MERNON TO FOR CE MALEVINE, (II) MISUSE, FAULT OR NEGLEGANCE OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WITHOUT THEY WERNON CONFORMINE AND TESTOR AND REVENT OF FORCE MALEVINE, (III) MISUSE, FAULT OR NEGLEGANCE OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WITHOUT THEY WERNON CONFORMINE PROPUNCTS, INARRES FOR WITHOUT THEY WERNON C

2