# Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight<sup>™</sup> 680

### **Product Details**

Size	500 µg	
Species Reactivity	Human	
Host/Isotype	Goat / IgG	
Class	Polyclonal	
Туре	Secondary Antibody	
Conjugate	DyLight™ 680	
Excitation/Emission Max	676/705 nm	
Immunogen	Human IgM	
Form	Liquid	
Concentration	0.5 mg/mL	
Purification	Antigen affinity chromatography	
Storage buffer	PBS, pH 6.8 to 7.4, with 0.2% BSA	
Contains	0.09% sodium azide	
Storage conditions	4° C	
RRID	AB_2556686	


Applications	Tested Dilution	Publications
Western Blot (WB)	1:5,000-1:20,000	-
Immunohistochemistry (IHC)	1:50-1:2,000	-
Immunocytochemistry (ICC/IF)	1:50-1:2,000	-
Flow Cytometry (Flow)	1:50-1:200	-
Immunoprecipitation (IP)	Assay-dependent	-

#### **Product Specific Information**

This antibody is cross-adsorbed and exhibits minimum reactivity to mouse and rat.


1

## Product Images For Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight™ 680

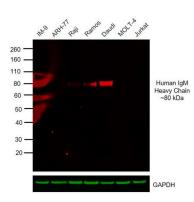


Human IqM Cross-Adsorbed Secondary Antibody (SA5-10106) in WB Western blot was performed using Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight 680 (Product # SA5-10106) and an ~80 kDa band corresponding to Human IgM Heavy Chain was observed in Raji, Ramos and Daudi but not in IM-9, ARH-77, MOLT-4 and Jurkat8203. Whole cell extracts (30 µg) of IM-9 (Lane 1), ARH-77 (Lane 2), Raji (Lane 3), Ramos (Lane 4), Daudi (Lane 5), MOLT-4 (Lane 6) and Jurkat (Lane 7) were electrophoresed using NuPAGE<sup>™</sup> 4-12% Bis-Tris Protein Gel (Product # NP0322BOX). Resolved proteins were then transferred onto a nitrocellulose membrane (Product # IB23001) by iBlot® 2 Dry Blotting System (Product # IB21001). The blot was probed with Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight 680 (Product # SA5-10106) (1:2000 dilution) and detected by fluorescence using the iBright FL1500 (Product # A44115). Raji, Ramos and Daudi are known to express IgM whereas IM-9 and ARH-77 express IgG and are negative for IgM. MOLT-4 and Jurkat, being T-cell lines, do not express immunoglobulins. (DOI: 10.1002/eji.1830100305; 10.3791/3573; 10.1016/0022-1759(94)00286-6; PMID: 566614).

Human IgM Cross-Adsorbed Secondary Antibody (SA5-10106) in ICC/IF

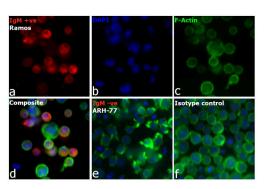

Immunofluorescence analysis of Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight680 was performed using log phase Ramos cells (IgM producing B-cell line). The cells were fixed with 4% paraformaldehyde for 10 minutes, permeabilized with 0.1% Triton<sup>™</sup> X-100 for 15 minutes and blocked with 2% BSA for 1 hour at room temperature. The cells were labeled with Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight680 (Product # SA5-10106) at 1:250dilution in 0.1% BSA, incubated at 4 degree celsius overnight (Panel a: red). Nuclei (Panel b: blue) were stained with Hoechst33342 (Product # H1399). F-actin (Panel c: green) was stained with Alexa Fluor™ 488 Phalloidin (Product # A12379,1:300 dilution). Panel d represents the merged image showing cytoplasmic (plasma membrane and golgi-body like) localization. Panel e represents ARH-77 (IgM non-producing B-cell line) which is a negative model for IgM expression.Panel f represents control cells with isotype control antibody to assess background. The images were captured at 40X magnification in Cell Insight CX7 LZR High-Content Screening (HCS) Platform (Product # CX7A1110LZR) (DOI:10.1002/eji.1830100305; 10.3791/3573; 10.1016/0022-1759(94)00286-6; PMID: 566614).

Human IgM Cross-Adsorbed Secondary Antibody (SA5-10106) in ICC/IF

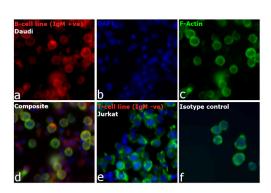

Immunofluorescence analysis of Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight680 was performed using log phase Daudi cells (Bcell line). The cells were fixed with 4% paraformaldehyde for 10 minutes, permeabilized with 0.1% Triton™ X-100 for 15 minutes and blocked with 2% BSA for 1 hour at room temperature. The cells were labeled with Goat anti-Human IgM Cross-Adsorbed Secondary Antibody, DyLight680 (Product # SA5-10106) at 1:250 dilution in 0.1%BSA, incubated at 4 degree celsius overnight (Panel a: red). Nuclei (Panel b: blue) were stained with Hoechst 33342 (Product# H1399). F-actin (Panel c: green) was stained with Alexa Fluor™ 488 Phalloidin (Product # A12379, 1:300 dilution). Panel d represents the merged image showing cytoplasmic (plasma membrane and golgi-body like) localization. Panel e represents Jurkat cells (T-cell line) which is a negative model for IgM expression. Panel f represents control cells with isotype control antibody to assess background. The images were captured at 40X magnification inCellInsightCX7 LZR High-Content Screening (HCS) Platform (Product # CX7A1110LZR) (DOI: 10.1002/eji.1830100305; 10.3791/3573; 10.1016/0022-1759(94)00286-6; PMID: 566614).

#### View more figures on thermofisher.com

2

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrate of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTIES, EXPRESS OR IMPLED, ARER CRANTED INCLUDING WITHOUT LIMITATION, IMPLED WARRANTIES OF MERCHANTABILITY. THRESS FOR ANY PARTICULAR PURPROSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS OUTING THE WARRANTY PERIOD IS UNITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AS THERE SON OBLIGATION TO FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE NON-CONFORMING PRODUCTS). A MANNER FOR EVENT OF FOR CEMALURE, (III) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE NON-CONFORMING PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unlass otherwise expressis stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo thrapeutic uses, or any type of consumption by or application to human or animals.